

General Information

Intended use:

FOB Turbilatex is a latex turbidimetric assay for **the quantitative detection of human haemoglobin (hHb) (faecal occult blood) in human stool samples.**

This assay is simple and widely applicable. Test results aid in a presumptive diagnosis of faecal occult blood (gastrointestinal bleeding).

For professional *in vitro* diagnostic use only.

FOB Turbilatex can be performed on every open chemistry analyser. Please follow the subsequent instructions in order to assure performance characteristics as describes in the instructions for use. This instruction has been validated by CerTest BIOTEC S.L Laboratories.

Additionally please read the "Instructions for use" for instructions on operating and programming user defined test.

Reagents:

Materials provided by CerTest BIOTEC S.L.:

Reagents	Quantity	Cat. reference
Turbidimetric reagents (R1 & R2) 200 Det/kit	R1: 2 vials, 2x22 mL. R2: 1 vial, 1x13 mL.	TL-022FB01 TL-022FB02
Auxiliary Reagents	Quantity	Cat. reference
Calibration kit	Calibrator: 6 vials, 6x1 mL.	TL-022FB70, TL-022FB71, TL-022FB72, TL-022FB73, TL-022FB74, TL-022FB75
Controls kit	Control C1, 2 vials, 2x1 mL/vial. Control C2, 2 vials, 2x1 mL/vial.	TL-022FB08 TL-022FB09
Sample diluent kit	4 vials, 4x125 mL/vial.	TL-022FB03E
Sample dilutions vials	1x2 mL/vial 1x2 mL/vial	MST-0005MF MST-0009F

Preparation of reagents:

R1 and R2 are ready to use.

Calibrators are ready to use.

Controls are ready to use.

Storage and stability

Kit components must be stored at temperature indicated on the label. Do not freeze.

Reagents are stable up to the expiration date printed on the label, always considering that reagent containers must be properly closed to avoid any contamination, must be kept away from the sunlight and conserved at temperature indicated on the label of each reagent.

Specimen:

Collect enough quantity of human stool samples. These samples should be collected in clean and dry containers (no preservatives or transport media). The samples can be stored completely dissolved in the refrigerator (2-8°C) for 3 days prior to testing. Homogenise stool samples as thoroughly as possible prior to preparation.

The sample dilution vial with diluted sample can be stored for 7 days in the refrigerator (2-8°C) prior to testing.

Use FOB Turbilatex stool collection tubes for sample collections described the instructions for use.

Assay procedure

Application parameter set up:

Specific analyzers settings for FOB Turbilatex must be programmed onto the analyzer, see below. For instructions, consult the TC220 (Tecom) analyzer manual and instructions for use provided with the kit.

Loading of reagents:

Load reagents according to the TC220 (Tecom) analyzer manual.

Calibration curve establishment:

A 6 point calibration curve can be established in TC220 (Tecom) analyzer. For instructions consult analyzer manual.

Calibration stability:

Calibrate the system at least once a month is extremely recommended. Recalibrate the system when reagent lot is change or when the controls are out of the assigned range given in the control labels and CoA.

QC controls:

FOB Turbilatex controls C1 and C2 must be assayed each day before running patient faecal sample extract to validate the calibration curve. The controls have assigned value ranges indicated on the label and certificate of analysis supplied. The control measurements must be within the indicated value range to obtain valid results for patient faecal extract. If the control values are out of range, follow next procedures: 1) Repeat QC control measurement, 2) Repeat calibration measurement.

Results:


The results are evaluated automatically by the analyzer and presented in ng hHb/mL.

Performance characteristics

The following results have been obtained during the validation of FOB Turbilatex on the TC220 (Tecom) analyzer.

Linearity:

FOB Turbilatex on TC220 (Tecom) analyzer using calibrator kit is linear in the calibration range of 0-1000 ng hHb/mL.

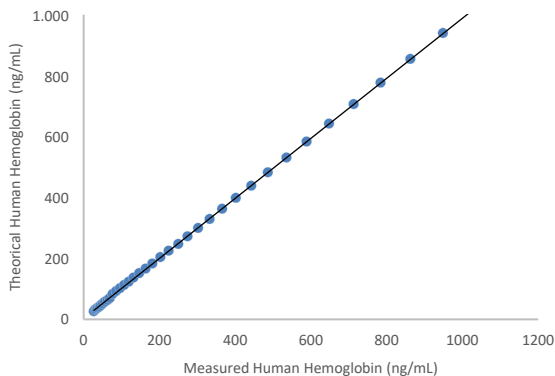

Measuring range:

FOB Turbilatex assay measuring range is 20-1000 ng hHb/mL on the TC220 (Tecom) analyser. Samples higher concentrated than 1000 ng hHb/mL must be diluted for proper quantification by the user, using additional sample buffer.

Prozone effect

Using the reported parameters, no prozone effect (hook effect) was observed up to 10 µg hHb/mL. Samples with Haemoglobin concentration of 16 µg/mL give a typical positive result 1000 ng hHb/mL.

Detection limit

Limit of detection (LOD): 10 ng hHb/mL. The lower limit of detection of FOB Turbilatex was determined on 20 samples and 2 sample replicates as the mean value + 2·SD.

Limit of quantification (LOQ): 20 ng hHb/mL. The lower limit of quantification is defined as the lowest actual amount of analysis that can be reliably detected; imprecision is < 20% as CV% on the TC220 (Tecom) analyzer.

Precision

FOB Turbilatex was tested with three different controls levels.

Method comparison

	Low (50 ng/mL)	Medium (250 ng/mL)	High (1000 ng/mL)
N	20	20	20
Mean (ng/mL)	49.9	232.8	947.2
SD (ng/mL)	2.3	9.8	21.1
CV (%)	5	4	2

Results obtained with FOB Turbilatex on the TC220 (Tecom) analyzer were compared with those obtained with EIKEN FOB Latex.

	Sensitivity	Specificity
FOB Turbilatex vs EIKEN FOB Latex	96%	>99%

Shipping damage











Please notify your distributor, if this product was received damaged.

NOTES

Please refer to the instruction for use for the detailed information about the test on the following:

Synthesis; Principle; Precautions; Reagents; Specimen collection and preparation; Interpretation of results and limitations.

Symbols key

	For <i>in vitro</i> diagnostic use only		Keep dry
	Consult instructions for use		Temperature limitation
	Catalogue number		Lot number
	Use by		Manufacturer
	Contains sufficient for <n> test	DIL	Sample diluent
	Keep out of the sunlight		

Manufacturer

CERTEST BIOTEC

Pol. Industrial Río Gállego II, Calle J, Nº 1, 50840,
 San Mateo de Gállego, Zaragoza (SPAIN)
 www.certest.es

TC220 (Tecom) /Application parameters

ASSAY PARAMETERS	
Std. No	6
R1	150 µL
Sample	30 µL
R2	20 µL
Others	NA
Reaction mode	Two points
Primary wavelength	510 nm
Secondary wavelength	None
Direction	Increase
Check time	300
Reaction time	10
Linear range	20-1000 ng/ml
CALIBRATION	
Calibration Method	Spline
Calibration set	6 calibrators
Blank	Calibrator 1 (0 ng/ml)
Calibrator 1	Calibrator 2 (50 ng/ml)
Calibrator 2	Calibrator 3 (100 ng/ml)
Calibrator 3	Calibrator 4 (250 ng/ml)
Calibrator 4	Calibrator 5 (500 ng/ml)
Calibrator 5	Calibrator 6 (1000 ng/ml)
STEPS	
Addition R1	0
Addition Sample	120
Incubation	
Addition R2	220
Blank Lecture	
Incubation	
Final lecture	



TC Series Automatic chemistry Analyzer

Test | Informes | QC | AJUSTES | Detener | Salir

Test Bioq. | Otros Parametros | Param. Calculados | Ajuste Perfiles | Secuencia Ensayo | Secuencia Impresion

Item: TP, ALB, TBA, GLU, CRE, UA, TC, TG, CK, LDH, AMY, Ca, Mg, Fe, NHS, MgO, CREO, P, PO, FTO, UREA, CLR, LACT, CALP2, **FOB**, PYL

Simbolo: FOB Nombre: FOB Cod. de barras: Param. Basicos Param. Calibracion

metodo: Two Points 1ª Long.0: 510 2ª long.0: NO Decimal: 1 Und: ng/ml
 (Rtvo.1+Rtvo.2)+S

Rtvo.1 Vol. (10-400ul): 150 Incub. Rtvo.: 120 Rtvo.1 Pos.: 4 Vol. Del envase: 20 ml 2ª Posicion: 0
 Rtvo.2 Vol. (0-400ul): 30 Rtvo.2 Pos.: 5 Vol. Del envase: 20 ml 2ª Posicion: 0

S Vol. (1.0-50.0ul): 20.0 Incub. Rtvo.: 100 Vel. de Mezcla: 0 Tiempos de reaccion: 10 Check Time(s): 300

Rango Linealidad: 20 — 1000 Correccion de resultados: $1.0000 X + 0.0000$ ratio de Dilucion (2-150): 10
 Absorbance Range: -0.5 — 3.5 Agotamiento de Sustrato: < 10.0 Vol. Muestra (3-70): 15
 m. de linealidad: 20 Vol. Total (120-500): 150
 o de valor Alto: Discrepancia = 10 Punto detectado en la Pendiente = 10 Ratio de Dilucion: 1:9
 test blanco Reactivo Blanco de Rtvo.: 0.0 — 10.0 Blanco Reactivo: 0.0000 todos ciclos dil.: 0
 todo el Vol.n de dilucion: 0

Fabricante: BIOLABO No. Lote: Caducidad: 30/06/2022

Añadir | Borrar | Guardar

MATERLAB Super User aaa 10/04/2019 11:38:20 Vers