

## General Information

### Intended use:

Transferrin Turbilatex is a latex turbidimetric assay for the quantitative detection of human transferrin (hTf) in human stool samples.

This assay is simple and widely applicable. Test results aid in a presumptive diagnosis of faecal occult blood (gastrointestinal bleeding).

For professional *in vitro* diagnostic use only.

Transferrin Turbilatex can be performed on every open chemistry analyser. Please follow the subsequent instructions in order to assure performance characteristics as describes in the instructions for use. This instruction has been validated by CerTest BIOTEC S.L Laboratories.

Additionally, please read the "Instructions for use" for instructions on operating and programming user defined test.

### Reagents:

#### Materials provided by CerTest BIOTEC:

Reagents	Quantity	Code
Turbidimetric reagents (R1 & R2) 200 Det/kit	R1: 2 vials, 2x22 mL R2: 1 vial, 1x13 mL	TL-022TF01 TL-022TF02
<b>Auxiliary Reagents</b>		
Calibration kit	Calibrator: 6 vials, 6x1 mL.	TL-022TF70, TL-022TF71 TL-022TF72 TL-022TF73 TL-022TF74 TL-022TF75
Controls kit	Control C1, 2 vials, 2x1 mL/vial. Control C2, 2 vials, 2x 1 mL/vial.	TL-022TF08 TL-022TF09
Sample diluent kit	4 vials, 4x125 mL/vial	TL-022UN03E
Sample dilutions vials	1x2 mL/vial 1x2.4 mL/vial	MST-0018MU MST-0019U

### Preparation of reagents:

R1 and R2 are ready to use.

Calibrators are ready to use.

Controls are ready to use

### Storage and stability

Kit components must be stored at temperature indicated on the label. Do not freeze.

Reagents are stable up to the expiration date printed on the label, always considering that reagent containers must be properly closed to avoid any contamination, must be kept away from the sunlight and conserved at temperature indicated on the label of each reagent.

### Specimen:

Collect enough quantity of human stool samples. These samples should be collected in clean and dry containers (no preservatives or transport media). The samples can be stored in the refrigerator (2-8°C) for 3 days prior to testing. For longer storage, maximum 6 months, the specimen must be kept frozen at -20°C. In this case, the sample will be totally thawed, and brought to room temperature (15-30°C) before testing. Freezing and thawing cycles are not recommended. Homogenise stool samples as thoroughly as possible prior to preparation.

The sample dilution vial with diluted sample can be stored for 7 days in the refrigerator (2-8°C) prior to testing.

Use Transferrin Turbilatex stool collection tubes for sample collections described the instructions for use.

## Assay procedure

### Application parameter set up:

Specific analyzers settings for Transferrin Turbilatex must be programmed onto the analyzer, see below. For instructions, consult the Alinity c-Series (Abbott) analyzer manual and instructions for use provided with the kit.

### Loading of reagents:

Load reagents according to the Alinity c Series (Abbott) analyzer manual.

### Calibration curve establishment:

A 6 point calibration curve can be established in Alinity c-Series (Abbott) analyzer. For instructions consult analyzer manual.

### Calibration stability:

Calibrate the system at least once a month is extremely recommended. Recalibrate the system when reagent lot is change or when the controls are out of the assigned range given in the control label and CoA.

### QC controls:

Transferrin Turbilatex controls C1 and C2 must be assayed each day before running patient fecal sample extract to validate the calibration curve. The controls have assigned value ranges indicated on the label and certificate of analysis supplied. The control measurements must be within the indicated value range to obtain valid results for patient fecal extract. If the control values are out of range, follow next procedures: 1) Repeat QC control measurement, 2) Repeat calibration measurement.

### Results:


The results are evaluated automatically by the analyzer and presented in ng hTf/mL.

## Performance characteristics

The following results have been obtained during the validation of Transferrin Turbilatex on the Alinity c-Series (Abbott) analyzer.

### Linearity:

Transferrin Turbilatex on Alinity c-Series (Abbott) analyzer using calibrator kit is linear in the calibration range of 0-250 ng hTf/mL.

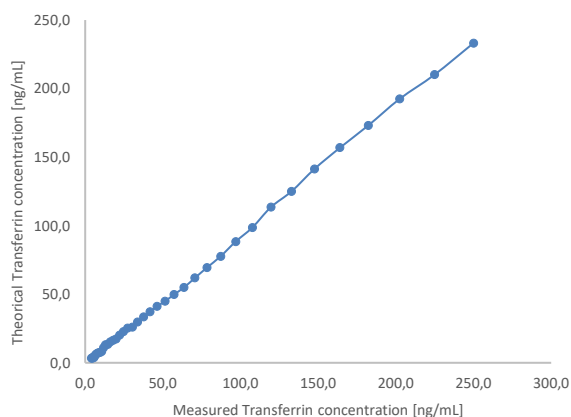

### Measuring range:

Transferrin Turbilatex assay measuring range is 2-250 ng hTf/mL on the Alinity c-Series (Abbott) analyser. Samples higher concentrated than 250 ng hTf/mL must be diluted for proper quantification by the user, using additional sample buffer.

### Prozone effect

Using the reported parameters, no hook effect was observed up to 10 µg hTf/mL. Samples with Transferrin concentration of 10 µg/mL give a typical positive result >250 ng hTf/mL.

### Detection limit

**Limit of detection (LOD): 1.4 ng hTf/mL. (\*)** The lower limit of detection of Transferrin Turbilatex was determined on 20 samples and 2 sample replicates as the mean value + 2 SD.

**Limit of quantification (LOQ): 2 ng hTf/mL (\*).** The lower limit of quantification is defined as the lowest actual amount of analysis that can be reliably detected; imprecision is < 20% as CV% on the Alinity c-Series (Abbott) analyzer.

**(\*) Data obtained by the analyser Biolis 24i (Tokyo Boeki)**

### Precision

Transferrin Turbilatex was tested with three different controls levels.

	Low (15 ng/mL)	Medium (80 ng/mL)	High (200 ng/mL)
<b>N</b>	20	20	20
<b>Mean (µg/g)</b>	14.8	83.1	202.4
<b>SD (µg/g)</b>	1.1	3.6	8.1
<b>CV (%)</b>	7.4	4.3	4.0

### Method comparison

Results obtained with Transferrin Turbilatex on the Alinity c-Series (Abbott) analyzer were compared with an immunochromatographic test (CerTest Transferrin, CerTest).

The results were as follows:

	Sensitivity	Specificity
<b>Transferrin Turbilatex vs CerTest Transferrin</b>	<b>95%</b>	<b>&gt;99%</b>

### Shipping damage

<b>IVD</b>	For in vitro diagnostic use only		Keep dry
	Consult instructions for use		Temperature limitation
<b>REF</b>	Catalogue number	<b>LOT</b>	Lot number
	Use by		Manufacturer
	Contains sufficient for <n> test	DIL	Sample diluent
	Keep out of the sunlight		

Please notify your distributor, if this product was received damaged.

### Symbols key

### Manufacturer

#### CERTEST BIOTEC

Pol. Industrial Río Gállego II, Calle J, Nº 1, 50840,  
San Mateo de Gállego, Zaragoza (SPAIN)  
www.certest.es

### NOTES

Please refer to the instruction for use for the detailed information about the test on the following:

**Synthesis; Principle; Precautions; Reagents; Specimen collection; Interpretation of results.**

**Alinity c-Series (Abbott) /Application parameters**

<b>ASSAY PARAMETERS</b>	
Std. No	6
R1	120 µL
Sample	10 µL
R2	30 µL
Others	NA
Reaction mode	End up
Primary wavelength	500 nm
Secondary wavelength	None
Direction	Increase
Self blank	21-21 cycle
Final lecture	34-36 cycle
Reaction time	10 min
Linear range	0-250 ng/ml
<b>CALIBRATION</b>	
Calibration Method	Linear
Calibration set	5 calibrators+Blank
Blank	Calibrator 1 (0 ng/ml)
Calibrator 1	Calibrator 2 (10 ng/ml)
Calibrator 2	Calibrator 3 (25 ng/ml)
Calibrator 3	Calibrator 4 (50 ng/ml)
Calibrator 4	Calibrator 5 (100 ng/ml)
Calibrator 5	Calibrator 6 (250 ng/ml)
<b>STEPS</b>	
Addition R1	
Addition Sample	
Incubation	
Addition R2	
Blank Lecture	Cycle 21
Incubation	
Final lecture	Cycle 34-36