

## General Information

### Intended use:

FOB Turbilatex is a latex turbidimetric assay for the quantitative detection of human haemoglobin (hHb) in human stool samples.

This assay is simple and widely applicable. Test results aid in a presumptive diagnosis of faecal occult blood (gastrointestinal bleeding).

For professional *in vitro* diagnostic use only.

FOB Turbilatex can be performed on every open chemistry analyser. Please follow the subsequent instructions in order to assure performance characteristics as describes in the instructions for use. This instruction has been validated by CerTest BIOTEC S.L Laboratories.

Additionally, please read the "Instructions for use" for instructions on operating and programming user defined test.

### Reagents:

#### Materials provided by CerTest BIOTEC S.L.:

Reagents	Quantity	Code
Turbidimetric reagents (R1 & R2) 200 Det/kit	R1: 2 vials, 2x15 mL R2: 1 vial, 1x6 mL	TL-022FB01C TL-022FB02C
<b>Auxiliary Reagents</b>		
Calibration kit	Calibrator: 6 vials, 6x1 mL.	TL-022FB70, TL-022FB71 TL-022FB72 TL-022FB73 TL-022FB74 TL-022FB75
Controls kit	Control C1, 2 vials, 2x1 mL/vial. Control C2, 2 vials, 2x 1 mL/vial.	TL-022FB08 TL-022FB09
Sample dilutions vials	1x2 mL/vial 1x2.4 mL/vial	MST-0018MU MST-0019U

### Preparation of reagents:

R1 and R2 are ready to use.

Calibrators are ready to use.

Controls are ready to use.

### Storage and stability

Kit components must be stored at temperature indicated on the label. Do not freeze.

Reagents are stable up to the expiration date printed on the label, always considering that reagent containers must be properly closed to avoid any contamination, must be kept away from the sunlight and conserved at temperature indicated on the label of each reagent.

### Specimen:

Collect enough quantity of human stool samples. These samples should be collected in clean and dry containers (no preservatives or transport media). The samples can be stored in the refrigerator (2-8°C) for 3 days prior to testing. Homogenise stool samples as thoroughly as possible prior to preparation.

The sample dilution vial with diluted sample can be stored for 7 days in the refrigerator (2-8°C) prior to testing.

Use FOB Turbilatex stool collection tubes for sample collections described the instructions for use.

### Assay procedure

#### Application parameter set up:

Specific analyzers settings for FOB Turbilatex must be programmed onto the analyzer, see below. For instructions, consult the ChemWell®-T (Awareness Technology Inc.) analyzer manual and instructions for use provided with the kits.

#### Loading of reagents:

Load reagents according to the ChemWell®-T (Awareness Technology Inc.) analyzer manual.

#### Calibration curve establishment:

A 6 points calibration curve can be established in ChemWell®-T (Awareness Technology Inc.) analyzer. For instructions consult analyzer manual.

#### Calibration stability:

Calibrate the system at least once a month is extremely recommended. Recalibrate the system when reagent lot is change or when the controls are out of the assigned range given in the control label and CoA.

#### QC controls:

FOB Turbilatex controls C1 and C2 must be assayed each day before running patient fecal sample extract to validate the calibration curve. The controls have assigned value ranges indicated on the label and certificate of analysis supplied. The control measurements must be within the indicated value range to obtain valid results for patient fecal extract. If the control values are out of range, follow next procedures: 1) Repeat QC control measurement, 2) Repeat calibration measurement.

#### Results:


The results are evaluated automatically by the analyzer and presented in ng hHb/mL.

## Performance characteristics

The following results have been obtained during the validation of FOB Turbilatex on the ChemWell®-T (Awareness Technology Inc.) analyzer using calibrator kit is linear in the calibration range of 0-1000 ng hHb/mL.

### Linearity:

FOB Turbilatex on ChemWell®-T (Awareness Technology Inc.) analyzer using calibrator kit is linear in the calibration range of 0-1000 ng hHb/mL.

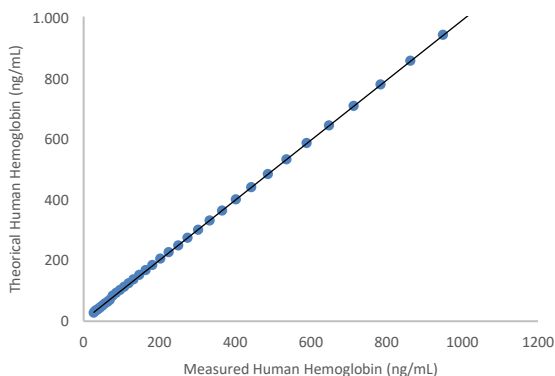

### Measuring range:

FOB Turbilatex assay measuring range is 15-1000 ng hHb/mL on the ChemWell®-T (Awareness Technology Inc.) analyser. Samples higher concentrated than 1000 ng hHb/mL must be diluted for proper quantification by the user, using additional sample buffer.

### Prozone effect

Using the reported parameters, no hook effect was observed up to 10 µg hHb/mL. Samples with Haemoglobin concentration of 10 µg hHb/mL give a typical positive result >1000 ng hHb/mL.

### Detection limit

**Limit of detection (LOD): 13 ng hHb/mL (\*).** The lower limit of detection of FOB Turbilatex was determined on 20 samples and 2 sample replicates as the mean value + 2·SD.

**Limit of quantification (LOQ): 15 ng hHb/mL (\*).** The lower limit of quantification is defined as the lowest actual amount of analysis that can be reliably detected; imprecision is < 20% as CV% on the ChemWell®-T (Awareness Technology Inc.) analyzer.

(\* Data obtained by the analyser Biolis 24i (Tokyo Boeki)

### Precision

FOB Turbilatex was tested with three different controls levels.

	Low (20 ng/mL)	Medium (80 ng/mL)	High (250 ng/mL)
<b>N</b>	20	20	20
<b>Mean (ng/g)</b>	22.0	84.2	254.2
<b>SD (ng/g)</b>	2.1	7.1	9.3
<b>CV (%)</b>	9	8	4

### Method comparison

Results obtained with FOB Turbilatex on the analyser Biolis 24i (Tokyo Boeki) analyzer were compared with those obtained with EIKEN FOB Latex.

	Sensitivity	Specificity
<b>FOB Turbilatex vs FOB Latex®</b>	<b>96%</b>	<b>&gt;99%</b>

### Shipping damage

Please notify your distributor, if this product was received damaged.

### Symbols key

	For in vitro diagnostic use only		Keep dry
	Consult instructions for use		Temperature limitation
	Catalogue number		Lot number
	Use by		Manufacturer
	Contains sufficient for <n> test	DIL	Sample diluent
	Keep out of the sunlight		

### Manufacturer

#### CERTEST BIOTEC

Pol. Industrial Río Gállego II, Calle J, N° 1, 50840,  
San Mateo de Gállego, Zaragoza (SPAIN)  
www.certest.es

### NOTES

Please refer to the instruction for use for the detailed information about the test on the following:

**Synthesis; Principle; Precautions; Reagents; Specimen collection; Interpretation of results.**

**ChemWell-T (Awareness) /Application parameters**


<b>ASSAY PARAMETERS</b>	
Std. No	6
R1	150 µL
Sample	20 µL
R2	30 µL
Others	NA
Reaction mode	Endpoint
Primary wavelength	505 nm
Secondary wavelength	None
Direction	Increase
Reagent blank lecture	10 seconds after R2 addition
Final lecture	300 seconds after blank lecture
Reaction time	10 min
Linear range	0-1000 ng/ml
<b>CALIBRATION</b>	
Calibration Method	Polinomial 4 <sup>th</sup> order
Calibration set	6 calibrators + Blank
Blank	Calibrator 1 (0 ng/ml)
Calibrator 1	Calibrator 2 (50 ng/ml)
Calibrator 2	Calibrator 3 (100 ng/ml)
Calibrator 3	Calibrator 4 (250 ng/ml)
Calibrator 4	Calibrator 5 (500 ng/ml)
Calibrator 5	Calibrator 6 (1000 ng/ml)
<b>STEPS</b>	
Addition R1	
Addition Sample	
Incubation	5 min
Addition R2	
Blank Lecture	
Incubation	5 min
Final lecture	

ChemWell-T Manager

Management Routines Lot# Registration QC Tracking Sample DB Settings Utilities View Help

Rack: EBDDD

Cooling: ON  
Current Temp: —

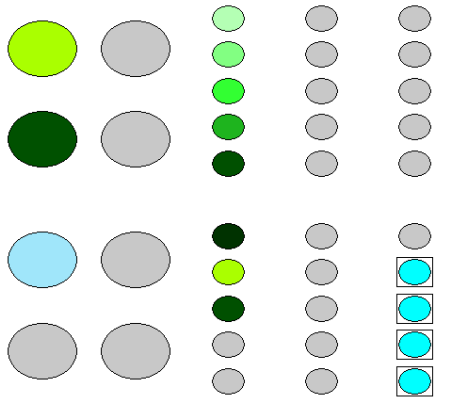

Clear Reagents Clear Samples

Load Home Save...

Cuvette Carrier

Heat: ON [37.0 °C]  
Current Temp: —

	1	2	3	4	5	
A	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
B	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
C	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
D	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
E	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
F	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
G	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear
H	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Clear

Verify Selected Clear All

Home

Prime Bottle Level: Unknown

Loading Instruction

Name/ID	Position	Volume
R1 FOB	Rack1:1	1600 ul
R2 FOB	Rack1:2	1200 ul
Cleaning Solution	Rack1:3	1500 ul
R-1	Rack1:38	+ 300 ul
R-2	Rack1:37	+ 300 ul
R-3	Rack1:36	+ 300 ul
R-4	Rack1:35	+ 300 ul

Cancel START

Layout Sample Calibration Test List Report

Ready NUM

Assay Editor for ChemWell-T

Assay View Substances Panels Indices Security Settings Help

Assay Name: FOB Turbilatex

Type: Chemistry Version: 1

Temperature: 37 °C

Assay Description:

Assay Mode: Polynomial 4th order  
Axis: Y = Abs, X = Conc

Assay Steps:

- Add Step: Add Reagent [R1 FOB], Volume [150.0 µL], One By One
- Add Step: Add Sample, Volume [20.0 µL]
- Edit Step: Clean Probe with [Cleaning Solution], Volume [100 µL]
- Incubate 00:05:00 [0h 5min 0sec]
- Add Step: Add Reagent [R2 FOB], Volume [40.0 µL], One By One
- Incubate 00:00:05 [0h 0min 5sec]
- Read with Primary Filter[505], Differential Filter[None]
- Incubate 00:05:00 [0h 5min 0sec]

Interpretation Settings:

Units: ng/mL

# of Decimals: 3

Special Groups

Normal Range:  Positive/Negative

Pos >= 50.000

Neg < 50.000

Reversed

Standards:

- St1: FOB Cal 0
- St2: FOB Cal 1
- St3: FOB Cal 2
- St4: FOB Cal 3
- St5: FOB Cal 4
- St6: FOB Cal 5

Curve Valid Time

Remove All

Blank:  Blank Used

Sample: Number of Replicates: 1

Controls:

- FOB Level 1
- FOB Level 2

For help, press F1

NUM